beurer

MG 10

ENGLISH

Vibrating Massage Ball

Instructions for use 2

ESPAÑOL

Pelota de masaje vibratoria

Manual de instrucciones..... 17

FRANÇAIS

Balle de massage vibrante

Mode d'emploi 34

**READ THIS MANUAL COMPLETELY AND CAREFULLY
BEFORE USING THIS PRODUCT**

Keep this manual in a safe location for future reference

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Questions or Comments? Call our US-based customer service toll free at 1-800-536-0366.

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ENGLISH

Table of contents

1. IMPORTANT SAFETY INSTRUCTIONS.....	3	6. Treatment Suggestions.....	9
2. Package Content.....	7	7. Care, Maintenance and Disposal	11
3. Parts and Controls.....	7	8. Technical specifications	12
4. Introduction	8	9. FCC Compliance Information.....	12
5. To Use	9	10. Warranty	14

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact customer service. See warranty for service contact.

Dear Valued Customer

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for application in the areas of heating, weight management, blood pressure, body temperature, pulse, gentle therapy, massage, and beauty.

Register your product today to receive further benefits:

- added protection for your purchase
- health and wellness tips
- new product information
- continual outstanding customer service and support.

Visit www.registerbeurer.com or contact customer service at 1-800-536-0366.

Please read these instructions for use carefully and keep them for later use. Make them accessible to other users and observe the information they contain.

With kind regards,
Your Beurer team

1. IMPORTANT SAFETY INSTRUCTIONS

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

READ THIS ENTIRE MANUAL, THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.

**WARNING**

WARNING indicates a hazardous situation which, if not avoided, could result in death or serious injury.

**CAUTION**

CAUTION indicates a hazardous situation which, if not avoided, may result in minor or moderate injury.

NOTICE

NOTICE addresses practices not related to personal injury, such as product and/or property damage

**WARNING – To reduce the risk of fire, electric shock, or serious personal injury:**

- Keep packaging material away from children. There is the risk of suffocation.
- Children must be supervised at all times while using the device.
- Never leave this device unattended, especially if children are nearby.
- Never use underneath a cover, such as a blanket, pillow, etc.
- Never use near gasoline or other highly flammable substances.
- Only use in dry interior rooms (e.g. never in the bath or sauna).
- Keep away from water or other liquids.
- Never reach for a device that is submerged in water.

- Only specialist personnel may perform repairs on electrical devices. Improper repairs may subject users to considerable danger. For repairs, please contact Customer Services or an authorized retailer.

CAUTION

- Before use, ensure that the device and accessories display no visible damage. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.
- Never open the housing.
- Keep the device away from pointed or sharp objects.
- If the device has been dropped or has suffered any other damage, it must no longer be used.
- The device must be switched off after every use and before cleaning.
- Do not use the massager:
 - On animals.
 - If you suffer from a medical abnormality or injury in the area to be massaged (e.g. herniated disc, open wounds).
 - Under any circumstances for massages in the area around the heart if you wear a pacemaker. If this is the case, consult your doctor regarding other parts of the body first.
 - On swollen, burnt, inflamed or damaged skin and body parts of any kind.
 - On cuts, capillaries, varicose veins, acne, couperose, herpes or other skin diseases.
 - During pregnancy.
 - On the face (eyes), on the larynx or other particularly sensitive body parts.
 - While sleeping.
 - After the consumption of medicines or alcohol (reduced cognition!).

- In a vehicle.
- In the public area.
- You should consult your doctor before using the massager, especially
 - If you suffer from a severe illness.
 - If you have a pacemaker, implant or other aids.
 - If you suffer from thrombosis.
 - If you suffer from diabetes.
 - If you are from suffering pain, the cause of which is unknown.

The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.

Battery Handling Safety Precautions

- Use only the size and type of batteries specified.
- Be sure to follow the correct polarity when installing the batteries. Reversed batteries may cause damage to the device.
- Do not mix different types of batteries together (e.g. Alkaline and Carbon-zinc or rechargeable batteries) or old batteries with fresh ones. Always replace batteries as a simultaneous set.
- If the batteries in the device are depleted or the device will not be used for a long period of time, remove the batteries to prevent damage or injury from possible battery leakage.
- Do not try to recharge batteries not intended to be recharged; they can overheat and rupture (follow battery manufacturer's directions.)
- Do not dispose of batteries in fire, batteries may explode or leak.
- Clean the battery contacts and also those of the device prior to battery installation.

- Remove discharged batteries from the product and dispose/recycle in compliance with all applicable laws.
- Keep batteries away from children and pets. Batteries may be harmful if swallowed. Should a child or pet swallow a battery, seek medical assistance immediately.

2. Package Content

1 x Vibrating Massage Ball MG 10

3 x 1,5V AAA Batteries

1 x Instruction Manual

3. Parts and Controls


On/Off button
(setting of the 2 intensity levels)

4. Introduction

This Vibrating Massage Ball performs targeted spot massage of tense muscle groups. It can also help to loosen connective tissue. The massage ball is comfortable to use because of its soft-touch surface. It features two intensity levels for both gentle and strong massage intensities.

Installing/Changing Batteries

1. Remove the packaging.
2. Check the device and batteries for damage.
3. Open the two halves of the massage ball by twisting them in opposite directions (twist the top half counterclockwise).
4. Insert three “AAA” batteries (included) into the battery compartment, ensuring the polarity (+/-) is correct.
5. Reassemble the two halves. When properly closed, the arrow markings on the two halves will align.

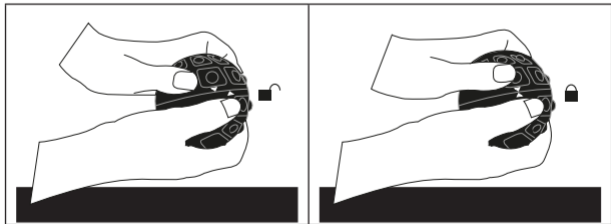5. To Use

CAUTION

The massage should always feel pleasant and relaxing. If you feel pain or discomfort, stop the massage, change your position, or decrease the contact pressure.

- You can switch on the device by pressing the On/Off button once. The massage ball starts to vibrate at the strongest intensity level.
- If you press the On/Off button again, the ball starts to vibrate at the weaker of the two intensity levels. Select one of the two intensity levels depending on your personal preferences.


Note

Treat each area (lower body, mid-body and upper body) three times a week for about 15 minutes each session. We recommend training all three areas, but you may choose to focus on one area. Perform each treatment slowly and carefully.

6. Treatment Suggestions


Calf

Sit on the floor and place one leg on the massage ball so that it is positioned on the bottom half of the calf. Support yourself with your hands while leaning back and lift your hips. Tense your stomach and bottom muscles and shift your lower body forwards and backwards. The calf muscles are stimulated.

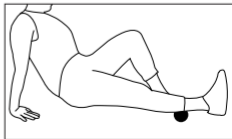
Thigh/bottom

Lie on your side on the floor and stretch both legs straight out. Bend one arm and support yourself on the floor with your elbow while keeping your other arm relaxed by your side. Place the massage ball under your hips. Now roll over the ball with slow movements. Vary the exercise by supporting yourself with both arms while leaning back; then roll your bottom over the ball with gentle pressing movements.

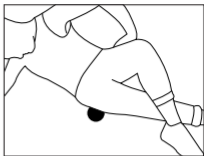Chest


Sit upright on the floor or on a chair and hold the ball in one hand. Now roll it over your chest in circular movements. In order to increase the pressure, you can also lean against a wall while standing up, with your front facing the wall.

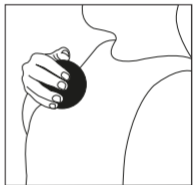Arm/shoulder

Lean against a wall while standing up, with your side facing the wall. Place the ball between your upper arm and the wall. Now apply gentle pressure on the ball and move up and down with your legs so that the ball rolls slowly along your arm. Increase the distance between your legs and the wall in order to increase the pressure on the ball.

Now slowly turn your back to the wall so that the ball moves along your shoulder area and the muscles in your shoulders are stimulated.

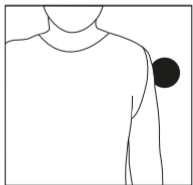
Back

Turn your back to the wall and place the ball on the right or left next to the spine between your back and the wall. Now move slowly up and down and roll the ball over your back.

CAUTION

Do not roll the ball directly over bones

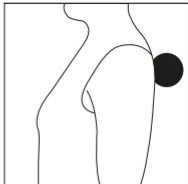7. Care, Maintenance and Disposal

NOTICE

- Switch off the device each time before cleaning.
- Clean the device using only the methods specified.
- Under no circumstances may water or liquid enter the device or accessories.
- Clean the device with a slightly damp cloth.
- Do not use abrasive cleaning products or stiff brushes.
- Ensure device is fully dry before use.
- Do not clean the device in a dishwasher.

Care

If you do not plan to use the device for a long period of time, we recommend that you store it in the original packaging in a dry environment, ensuring it is not weighed down by additional items.

Disposal

Observe the local regulations for material disposal. Dispose of the device in accordance with local regulations. If you have any questions, please contact the local authorities responsible for waste disposal.

8. Technical specifications

Type	MG 10
Weight	Approx. 5,7 oz (162 g)
Battery	3 x 1.5 V  AAA (LR03) battery
Device dimensions	Diameter approx. 3.15 in (8 cm)

9. FCC Compliance Information

Vibrating Massage Ball MG 10

Responsible Party – U.S. Contact Information:

Beurer North America LP

1 Oakwood Boulevard, Suite 255

Hollywood, FL 33020

United States

1-800-536-0366

info@beurer.com

FCC Compliance Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) this device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

10. Warranty

Limited Lifetime Warranty For Original Purchaser

Your Beurer Vibrating Massage Ball, model MG 10, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer Vibrating Massage Ball, model MG 10, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at 1-800-536-0366 or at info@beurer.com to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn batteries; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the

period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer. This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state. For more information regarding our product line in the USA, please visit: www.beurer.com

Distributed by:
Beurer North America LP
1 Oakwood Boulevard, Suite 255
Hollywood, FL 33020, USA
www.beurer.com

Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.

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